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English

When using this unit for the first time, please read the instructions and guidelines below carefully:


Power supply: 2 AA batteries Power: 2W Net weight range: 375g

Outside dimension: W130 x D900 x H85 (mm)

The battery is not included in the product.

Product Components:

1. button
2. Battery Cover
3. Battery box
4. control
5. Battery Bag
6. jack

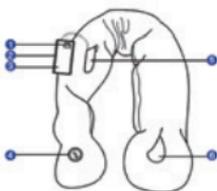Suggested use: neck, shoulder, leg, arm, waist.

Set the switch to "0" and open the bag to disconnect the power supply from the battery box.

No drag the wire during use or may occur a power interruption.

pescoço

ombro

braço

perna


cintura

• Cleaning method:

Clean the product with a damp cloth, only when it is dirty

(Warning: Do not wash in a washing machine)

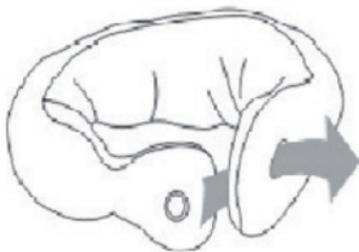The most common defects and how to solve them

Defect	Possible Cause	How to solve it
It does not work when connected	Incorrect positioning of the batteries	Replacing batteries
It does not work when connected	The power supply is not connected	Connect the power supply
Unusual noise when connected	Problem in engine vibration	Send an authorized repair
For working unexpectedly	Disconnected or damaged wiring	Send an authorized repair


How to solve a problem during an emergency:

- An unauthorized person not allowed to disassemble, repair or modify the product if defective.
- Call for service or contact your nearest distributor if you have any problems with the product


Use:

Open the bag and slide the battery cover on the side of the battery box (see image at right), place two AA batteries into position and slide the cover back.

Passes to the B side of the box and place the switch in the "ON" position, then place the battery box back in the bag and close the zipper. You can turn on the track taken when using it and fixed it at the massage you want (such as neck, shoulder, leg, arm or waist). Select the method of massage you want to control.

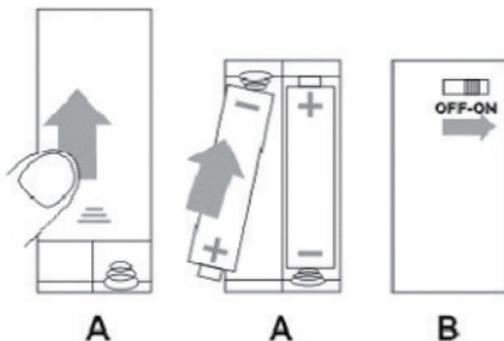- There are two intensities to select.
- Turn on switch, “I” means low speed, vibration is lower and slower. “II” means higher speed, vibration is stronger and faster. “O” means it is connected (see box at right). You can select different speeds according to their will.
- It is designed to fully cover the part that are massaged: you can simply plug the strip into the outlet and adjust the cable to the appropriate length to cover the neck, shoulder, leg, arm or waist completely. The internal vibrator will relax the muscles. The wireless design will offer you the comfort of getting to where you want.

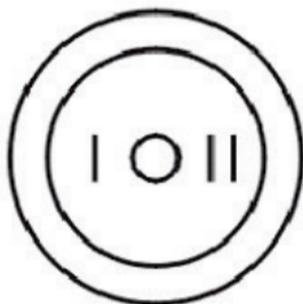- Usage time suggested: it is preferable to be 30 minutes. Take a break of 15 if you want to use it again.

www.multilaser.com.br

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MULTILASER

Massageador Wrap Massagem